




粒子径分布測定結果

Microtrac Version 11.0.0-246W

光学台	MT3300(LOW-WET)	サンプルID1	Face powder
レコードNo.	9	サンプルID2	
測定日付	2025/04/04	タイトル	YTC-D-2-16-2
測定時刻	10:08	シリアルNo.	S6453 (3328)
ファイル名	C:\Users\¥Microtrac¥Documents¥NMA¥Munekawa¥Munekawa.dms2		

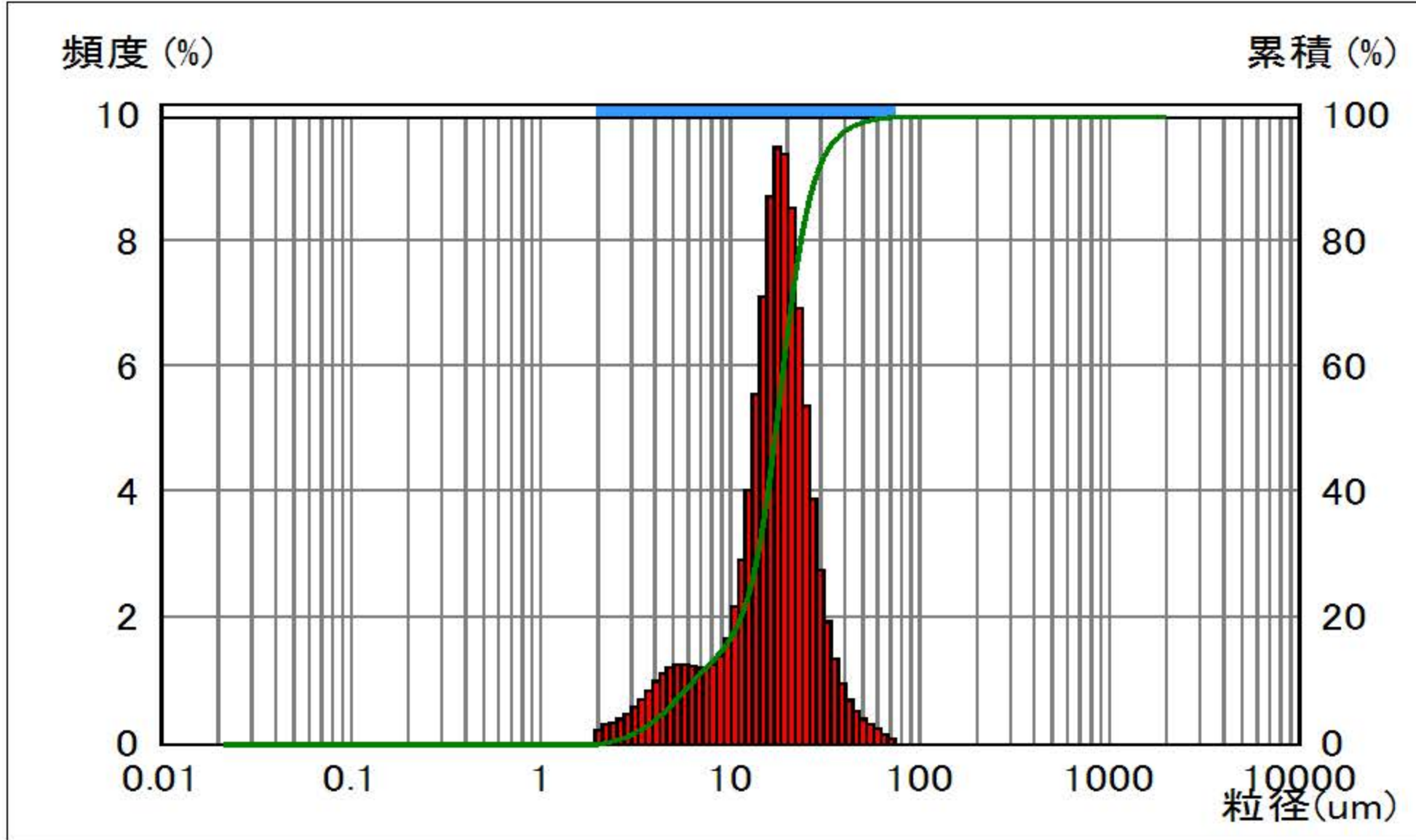

要約データ	値
MV(um)	17.888
MN(um)	4.4215
MA(um)	12.639
CS	0.47471
SD(um)	7.6388

ピーク(um)	頻度(%)	幅(um)
17.368	100.00	15.278

CH	粒径(um)	頻度(%)	累積(%)	CH	粒径(um)	頻度(%)	累積(%)	CH	粒径(um)	頻度(%)	累積(%)	CH	粒径(um)	頻度(%)	累積(%)
1	2000	0.00	100.00	34	114.1	0.00	100.00	67	6.541	1.21	10.47	100	0.375	0.00	0.00
2	1826	0.00	100.00	35	104.7	0.00	100.00	68	5.998	1.23	9.26	101	0.344	0.00	0.00
3	1674	0.00	100.00	36	95.96	0.00	100.00	69	5.500	1.23	8.03	102	0.315	0.00	0.00
4	1535	0.00	100.00	37	88.00	0.00	100.00	70	5.044	1.18	6.80	103	0.289	0.00	0.00
5	1408	0.00	100.00	38	80.70	0.00	100.00	71	4.625	1.08	5.62	104	0.265	0.00	0.00
6	1291	0.00	100.00	39	74.00	0.05	100.00	72	4.241	0.96	4.54	105	0.243	0.00	0.00
7	1184	0.00	100.00	40	67.86	0.11	99.95	73	3.889	0.81	3.58	106	0.223	0.00	0.00
8	1086	0.00	100.00	41	62.23	0.20	99.84	74	3.566	0.67	2.77	107	0.204	0.00	0.00
9	995.6	0.00	100.00	42	57.06	0.28	99.64	75	3.270	0.55	2.10	108	0.187	0.00	0.00
10	913.0	0.00	100.00	43	52.33	0.36	99.36	76	2.999	0.44	1.55	109	0.172	0.00	0.00
11	837.2	0.00	100.00	44	47.98	0.48	99.00	77	2.750	0.36	1.11	110	0.158	0.00	0.00
12	767.7	0.00	100.00	45	44.00	0.66	98.52	78	2.522	0.30	0.75	111	0.145	0.00	0.00
13	704.0	0.00	100.00	46	40.35	0.93	97.86	79	2.312	0.27	0.45	112	0.133	0.00	0.00
14	645.6	0.00	100.00	47	37.00	1.32	96.93	80	2.121	0.18	0.18	113	0.122	0.00	0.00
15	592.0	0.00	100.00	48	33.93	1.91	95.61	81	1.945	0.00	0.00	114	0.111	0.00	0.00
16	542.9	0.00	100.00	49	31.11	2.74	93.70	82	1.783	0.00	0.00	115	0.102	0.00	0.00
17	497.8	0.00	100.00	50	28.53	3.87	90.96	83	1.635	0.00	0.00	116	0.094	0.00	0.00
18	456.5	0.00	100.00	51	26.16	5.36	87.09	84	1.499	0.00	0.00	117	0.086	0.00	0.00
19	418.6	0.00	100.00	52	23.99	6.90	81.73	85	1.375	0.00	0.00	118	0.079	0.00	0.00
20	383.9	0.00	100.00	53	22.00	8.51	74.83	86	1.261	0.00	0.00	119	0.072	0.00	0.00
21	352.0	0.00	100.00	54	20.17	9.37	66.32	87	1.156	0.00	0.00	120	0.066	0.00	0.00
22	322.8	0.00	100.00	55	18.50	9.47	56.95	88	1.060	0.00	0.00	121	0.061	0.00	0.00
23	296.0	0.00	100.00	56	16.96	8.70	47.48	89	0.972	0.00	0.00	122	0.056	0.00	0.00
24	271.4	0.00	100.00	57	15.56	7.10	38.78	90	0.892	0.00	0.00	123	0.051	0.00	0.00
25	248.9	0.00	100.00	58	14.27	5.54	31.68	91	0.818	0.00	0.00	124	0.047	0.00	0.00
26	228.2	0.00	100.00	59	13.08	4.01	26.14	92	0.750	0.00	0.00	125	0.043	0.00	0.00
27	209.3	0.00	100.00	60	12.00	2.90	22.13	93	0.688	0.00	0.00	126	0.039	0.00	0.00
28	191.9	0.00	100.00	61	11.00	2.15	19.23	94	0.630	0.00	0.00	127	0.036	0.00	0.00
29	176.0	0.00	100.00	62	10.09	1.65	17.08	95	0.578	0.00	0.00	128	0.033	0.00	0.00
30	161.4	0.00	100.00	63	9.250	1.38	15.43	96	0.530	0.00	0.00	129	0.030	0.00	0.00
31	148.0	0.00	100.00	64	8.482	1.23	14.05	97	0.486	0.00	0.00	130	0.028	0.00	0.00
32	135.7	0.00	100.00	65	7.778	1.18	12.82	98	0.446	0.00	0.00	131	0.026	0.00	0.00
33	124.5	0.00	100.00	66	7.133	1.17	11.64	99	0.409	0.00	0.00	132	0.023	0.00	0.00

累積(%)	粒径(um)
10.00	6.3305
20.00	11.284
30.00	13.925
40.00	15.762
50.00	17.368
60.00	19.029
70.00	20.921
80.00	23.448
90.00	27.874
95.00	32.920

測定回数	Avg/2	測定上限(um)	2000	粒子名	SiO2	溶媒名	WATER
測定時間(秒)	30	測定下限(um)	0.021	粒子透過性	透過	溶媒屈折率	1.333
分布表示	体積	残分比(%)	0.00	粒子屈折率	1.48	フィルタ	標準
粒径区分	標準	通過分比(%)	0.00	粒子形状	非球形	DV値	0.0980
計算モード	MT3000II	残分比設定	無効			透過率(TR)	0.924
						拡張フィルタ	無効